mamajoo®

EN Mamajoo Electronic Single Breast Pump

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www.mamajoo.com

Dear Consumer,

Thank you for choosing Mamajoo products. Mamajoo products meet the stringent EU standards as they are designed and produced with highest quality BPA free raw materials combined with latest high technology, by considering health and safety. For more information and your questions, you can find our contact details at the end of this instructions leaflet. Please keep this leaflet for future references.

USER INSTRUCTIONS / WARRANTY**ATTENTION**

1. Please read all instructions carefully before use.
2. Do not clean or sterilize the motor unit, its components, electrical adapter, air tubing and their inlets because they are not in contact with breast milk. Before using for the first time, remove all parts of your pump that come in contact with your breast milk, wash them thoroughly with warm soapy water, rinse well and sterilize using Mamajoo sterilizers or by boiling them in water for 5 minutes. Clean all parts after each use and sterilize all parts before each use. To prevent damage to the pump components, do not sterilize longer than recommended.
3. NEVER expose parts that do not come into contact with breast milk, such as the motor unit and its components, electrical adapter, air tubing and inlets to water or moisture, do not clean with a wet cloth.
4. Avoid cleaning with anti-bacterial detergents and cleaners that may cause damage or cracking on the product. In such a case, please do not use damaged pieces.
5. Use Mamajoo Electronic Single Breast Pump with ONLY the original adapter. Check if it turns off correctly after each use and unplug the power cord. Thus, you save energy and extend the life of the product.
6. Do not exceed the capacity of the breast milk container to avoid motor unit damages. If you get more than 120 ml breast milk when you pump, please use a larger feeding bottle.
7. Keep out of reach of children and babies.
8. Extra hard and fast milk expression may cause pain so please ensure the level is right for you. We do not recommend you to use your breast pump at the last suction level except in special cases such as Mastitis.
9. Please check before each use that the comfort cushion fits correctly to the funnel, otherwise milk can leak out.
10. When a problem occurs, unplug the power cord. In case of any requirements of service, please contact us to direct you to our Authorized Technical Service.
11. There should be strict compliance with the instructions. Our company will not be held responsible for any damage caused by improper use.
12. Do not use the breast pump in the bath or in the shower.

13. Do not place or store the breast pump near a tub or sink. Do not reach for a product that has fallen into the water, unplug immediately.
14. Always unplug the breast pump immediately after use.
15. This product is not intended for use by people (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning the use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance
16. Do not allow children or pets to play with the motor unit, the adapter or product parts as these parts may cause choking danger.
17. Inspect the breast pump, including the adapter, for signs of damage before each use. Do not use the breast pump when the adapter or plug is damaged, does not work properly or has been dropped into the water.
18. There are no user serviceable parts inside the motor unit of the breast pump. Do not attempt to open, service or repair the breast pump yourself. If you do so, your guarantee becomes invalid.
19. Keep the adapter and the silicone tubes away from heated surfaces. Only use the adapter indicated on the appliance.
20. Do not attempt to remove the pump body from your breast while under vacuum. Switch off the appliance and break the seal between your breast and the pump funnel with your finger to remove the pump from your breast.
21. Always start milk expression at massage mode first level to stimulate let-down reflex. Do not continue pumping for more than 5 minutes at a time if you do not succeed in expressing any milk. Try to express at another time during the day. If the process becomes very uncomfortable or painful, stop using the pump and consult your breastfeeding advisor.
22. Check the product before each use and replace any worn parts. Expressing milk with worn out parts may result motor unit damage due to milk leakage into the motor.
23. If the valve is not cleaned well, expressing milk properly will not be possible. After each use, remove the valve and clean it well from any residual milk.

MAINTENANCE

1. Clean the motor unit ONLY with a dry cloth. Do not use abrasive or anti-bacterial detergents.
2. Check regularly that the air tube connections are made correctly and the valve is not damaged.
3. The pump body and the breast milk containers can be washed in the dishwasher top rack. However a color change can occur due to food colours. So, we recommend you to wash them separately.

The air tube CAN NOT be washed or sterilized. The water left in the air tube might destroy the motor unit.

TECHNICAL FEATURES

- The advanced intelligent microprocessor technology provides automatic control.
- Processes and signs can be easily monitored from LED indicators.


Product Code: MMJ4401

Adaptor: Input: 100- 240V ~ 50/60Hz Output: 5-6V 1A DC

Run Time: 0-30 minutes

Weight: 0.150 kg

Electronic Single Breast Pump product components

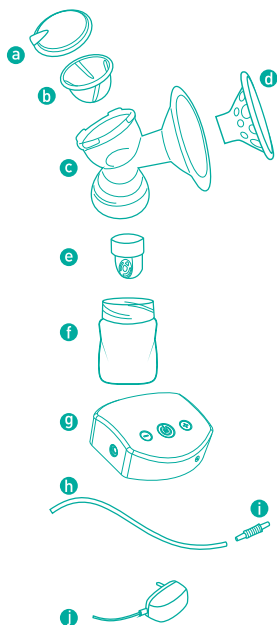
- a) Cylinder cover
- b) Silicone diaphragm
- c) Pump body
- d) Silicone comfort cushion
- e) Valve
- f) Breast milk container (160 ml)
- g) Motor unit
- h) Air tube
- i) Air tube – Motor unit insert
- j) AC/DC adapter

**NEVER sterilize these parts:
g, h, i, and j**

HOW TO ASSEMBLE AND LED INDICATOR MANUAL

Before starting to assemble, please wash your hands.

Before using wash and sterilize the mechanical pieces of the breast pump.

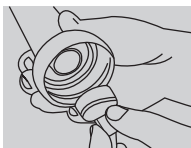
Push the valve to the pump body tightly.

Place the silicone diaphragm onto the pump body.


Place the cylinder cover onto the pump body.

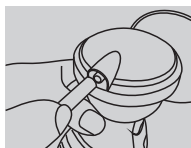
Insert the end of the air tube into the air port of the cylinder cover.

Screw breast milk container to the pump body.

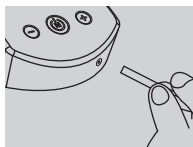Insert the air tube into the port of the motor unit

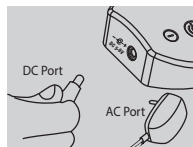Insert the metal plug at the end of adapter into the adapter input on the motor unit (5V-6V) and plug in the adapter (110V-220V).

Mamajoo Electronic Breast Pump

2 Different Message Mode Options

5 Different Suction Mode Options

Control Panel

- a) Massage Mode
- b) Suction Mode
- c) On/Off and Massage/Suction Mode Button
- d) Level Decrease Button
- e) Level Increase Button
- f) LED Level Indicator

OPERATION

1. Make sure the breast pump is assembled correctly.
2. Sit on a comfortable chair and relax, hold the pump with your hand which you feel comfortable with. Do not use breast pump in leaning position.
3. Center the pump funnel with silicone comfort cushion over your nipple. Make sure it seals your breast completely.
4. Plug in and press 'On/Off' button to start suction. We recommend you to start with massage mode to stimulate your milk flow.
5. Mamajoo Breast Pumps mimic baby's natural suckling rhythm. When you press the On/Off button of Mamajoo Electronic Single Breast Pump one time, LED indicator at the upper left lights up for the massage mode, if you press the button for a second time LED indicator at the upper right lights up for the suction mode. The massage mode is to stimulate the "let-down" reflex that makes the milk come in, as the baby does when starts breastfeeding. At the first 2 levels of the massage mode; you will hear a long vacuum shot after each 10 short shots and when your milk starts to flow, press the "On/Off" button again to switch to suction mode. You can select the desired expression power between 3-7 levels and with "+" or "-" buttons you can adjust suction power. Mothers with nipple problems and sensitivities can also express milk at massage mode. When the expression is finished, press "On/Off" button for 2 seconds and unplug it. Take out the breast milk container from the pump.
6. Mamajoo Electronic Breast Pump automatically shuts off after it runs 30 minutes, thanks to its intelligent micro-processor technology.

WARNING

The product should be moved by considering the instructions on its box, by protecting from high temperatures, fire and impacts in its original package during transportation and shipping.

WARRANTY

1. This product is under guarantee of our company for 2 years after date of purchase for the first user. Please contact us directly or the retailer you purchased the product for warranty claims.
2. This warranty does not cover wear & tear results occur from regular use, defects or damages caused by misuse, abuse and improper storage, failure to comply the manual instructions, malfunctions from using improper components.
3. The responsibility of the usage the product belongs to the user.
4. Please keep your receipt as a proof of purchase for any warranty claim.

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Tarih - İmza - Kaşe :

AWARDS / HONOURS

Mamajoo Electronic USB Breast Pumps won the **“Mothers 1st Choice”** award in 2017

Mamajoo feeding bottles are used by Lufthansa to provide their passengers since 2016.

Mamajoo, won the **“Limitless Export Champs”** award in 2016 Export Summit.

Mamajoo 5 in 1 Steam Steriliser & Dryer is chosen **“The Best Product”** in its category by German experts (www.expertentesten.de) since 2016.

Mamajoo 3 in 1 Steam Steriliser & Heater has been placed among the best products by German experts (www.expertentesten.de) in 2017.

Mamajoo feeding bottles were chosen as **“The Best Feeding Product”** in a survey conducted by the Baby&Kids Store Magazine in 2015.

Mamajoo feeding bottles were given the innovative **“Product & Design”** award in 2012 by PAGEV.

According to a search by Marketing Turkey, Mamajoo is awarded as **“The Most Prestigious Brand of The Year”** for 2019 in the Baby Care and Nutrition Category.

Mamajoo brand and manufacturer group company Mamajoo Health Products became the first company to receive the TSE Covid-19 Safe Production Certificate in its field.

Mamajoo Health Products have been awarded the HALAL Certificate approved by DSR Certification in the category of mother and baby products (Nutrition, breastfeeding and Care products) in 2020

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mamajoo GmbH